

CytoGrow™ Recombinant IL-17A (Human)


Catalog Numbers: CG121-A, CG121-B, CG121-C

Revision 1.0

Specifications

Synonym	
Species	Human
Size	10µg/50µg/1mg
Tag	Tag Free
Purity	≥ 95%
Endotoxin	≤10 EU/mg by the LAL method
Expression System	CHO
Expression Region	Ile20-Ala155
pH	7.0-8.0
Appearance	White Powder or clear liquid
Formulation	Lyophilized from a 0.22 µm-filtered solution containing PBS, 5% mannitol and 0.01% Tween 80, pH 7.4
Shipping	US - Ships overnight
Reconstitution	Recommended to redissolve in sterile deionized water.
Storage and Stability	<p>Lyophilized state: 36 months at -20°C to -80°C</p> <p>After reconstitution and under sterile conditions. 6 months at -20°C to -80°C 7-10 days at 2°C to 8°C</p> <p>Store in a manual defrost freezer and avoid repeated freeze-thaw cycles.</p>

Bioactivity

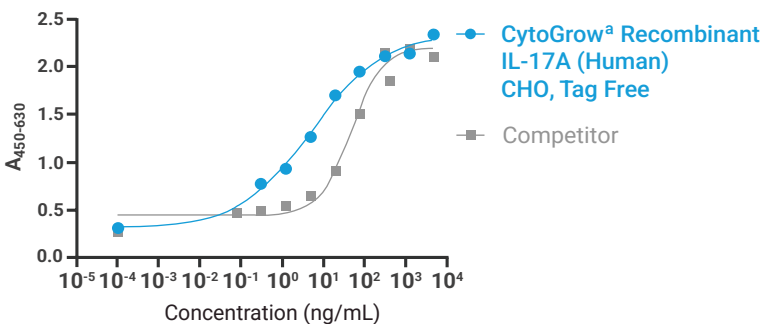

Measured by its ability to induce IL-6 secretion by NIH-3T3 mouse embryonic fibroblast cells. The expected ED₅₀ is less than or equal to 0.5 ng/mL.

Reconstitution

1. Reconstitute protein to concentrations of 100 µg/mL or higher.
2. Use sterile deionized water to solubilize the lyophilized protein.
3. Gently tap glass vial to help collect the powder to the bottom.
4. Add desired amount of buffer, observe after 2 min. Gentle agitation by hand may be required.
5. Prepare aliquots of solution. Freeze at -20°C to -80°C to avoid repeated freeze thaw cycles.